

# The new KleenSpec® Disposable Vaginal Specula and Cordless Illumination System

## Ordering Information

### KleenSpec® Vaginal Specula—Premium 590 Series

- 59000 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Small (24/box; 4 boxes/case)
- 59001 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Medium (24/box; 4 boxes/case)
- 59004 KleenSpec® 590 Series Premium Disposable Vaginal Specula, Large (19/box; 4 boxes/case)

### KleenSpec® Box Series

- 59000-B KleenSpec® 590 Series Premium Disposable Vaginal Specula, Small (24/box)
- 59001-B KleenSpec® 590 Series Premium Disposable Vaginal Specula, Medium (24/box)
- 59004-B KleenSpec® 590 Series Premium Disposable Vaginal Specula, Large (19/box)

### KleenSpec® Vaginal Specula—Premium 590 Series w/Smoke Tube

- 59005 KleenSpec® 590 Series Premium Disposable Vaginal Specula with Smoke Tube, Small (12/box; 4 boxes/case)
- 59006 KleenSpec® 590 Series Premium Disposable Vaginal Specula with Smoke Tube, Medium (12/box; 4 boxes/case)

### KleenSpec® Disposable Vaginal Specula Illumination System

- 79910 Complete KleenSpec® Disposable Vaginal Specula Cordless Illuminator with Charging Station
- 79915 Same as above, for export purposes—specify country and voltage; conforms to IEC 601-1
- 79920 Complete KleenSpec® Disposable Vaginal Specula Cordless Illuminator with Charging Station and Auxiliary Power Cord
- 79925 Same as above, for export purposes—specify country and voltage; conforms to IEC 601-1
- 79900 KleenSpec® Disposable Vaginal Specula Cordless Illuminator
- 73910 Charging Station for Cordless Illuminator;  
For export purposes: 73912 (European), 73914 (UK), 73916 (Australia)
- 75910 Auxiliary Power Cord Accessory for Cordless Illuminator;  
For export purposes: 75912 (European), 75914 (UK), 75916 (Australia)

In addition, we offer a complete line of Vaginal Specula and Accessories to meet your needs:

### KleenSpec® Vaginal Specula—Premium 580 Series w/Sheath

- 58000S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Small (20/box; 4 boxes/case)
- 58001S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Medium (20/box; 4 boxes/case)
- 58004S KleenSpec® 580 Series Premium Disposable Vaginal Specula with Sheath, Large (15/box; 4 boxes/case)

### KleenSpec® Box Series

- 58000S-B KleenSpec® 580 Series Premium Disposable Vaginal Specula w/Sheath, Small (20/box)
- 58001S-B KleenSpec® 580 Series Premium Disposable Vaginal Specula w/Sheath, Medium (20/box)

### KleenSpec® Vaginal Specula—Economy Series

- 58600 KleenSpec® 586 Series Economy Disposable Vaginal Specula, Small (25/box; 4 boxes/case)
- 58601 KleenSpec® 586 Series Economy Disposable Vaginal Specula, Medium (25/box; 4 boxes/case)
- 78600 Expandable Illuminator (6/box; 6 boxes/case)—Can be sterilized 50 times before disposal


### KleenSpec® Disposable Vaginal Specula Illumination System

- 78810 KleenSpec® Illumination System—4.7 V Halogen Light Handle with on/off switch and transformer. 110-130 V, 60 Hz for North America. UL Listed, CSA Certified
- 78815 Same as above, for export purposes—Specify country and voltage; conforms to IEC 601-1
- 78800 Permanent 4.7 V Halogen Light Handle with On/Off Switch
- 73410 Transformer with 5' Cord—110-130 V, 50 Hz. UL Listed, CSA Certified
- 73315 Same as 73410, but for export purposes—Specify country and voltage; conforms to IEC 601-1
- 08800-U 4.6 V to 4.7 V Halogen Replacement Lamp for 78810 Illumination System
- 07800-U 6 V Halogen Replacement Lamp for 78010 Illumination System

### KleenSpec® Storage Systems

- 78820 KleenSpec® Illumination Storage System
- 58014 KleenSpec® Vaginal Speculum Dispenser

Patents Pending

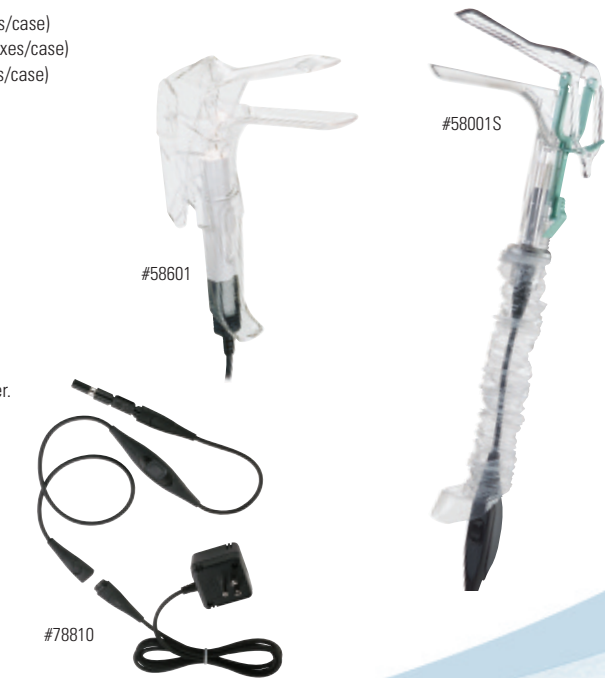

4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA  
(p) 800.535.6663 (f) 315.685.3361 www.welchallyn.com

WelchAllyn®

Advancing Frontline Care™

© 2007 Welch Allyn SM2582 Rev E



## Introducing the new KleenSpec® Disposable Vaginal Specula and Cordless Illumination System

An advancement in Women's Health technology that's nothing short of brilliant.



WelchAllyn®

Advancing Frontline Care™



# Now you can have better visualization of the exam area without the hassle of cords.

The new KleenSpec® Disposable Vaginal Speculum with the innovative Cordless Illumination System offers greater convenience for you, added comfort for your patients.

The Welch Allyn KleenSpec® Vaginal Specula and Illumination System has been totally redesigned to provide the key features you've come to expect—plus a number of important and exciting enhancements.

With so many new time-saving features and benefits, the new KleenSpec® Disposable Vaginal Specula and Cordless Illumination System are a must for every office and hospital.



## The new 590 Series KleenSpec® Vaginal Speculum has all of the quality features you rely upon and more.

- > **Maneuverability**—new wider handle provides improved ergonomics, better balance, and it is much easier to manipulate during exams.
- > **Superior Illumination**—newly designed light pipe provides a homogeneous spot and eliminates glare back—all of the light is projected forward right where you need it.
- > **Comfortable**—smooth, molded acrylic will not pinch or bind tissue for greater patient comfort during the exam.
- > **Disposable**—Welch Allyn's single-use speculum eliminates the time and expense of reprocessing and reduces the chance of cross-contamination.
- > **Safe and Dependable**—stronger grip and the yoke is more flexible than low-cost substitutes.
- > **Compatibility**—the new 590 Series KleenSpec® Vaginal Speculum is compatible with the current illumination system (#78810) as well as the new illumination system (#79910).



590 Series KleenSpec® Vaginal Speculum

Available in small (#59000), medium (#59001) and large (#59004) specula

## The new Welch Allyn KleenSpec® Disposable Vaginal Specula 790 Series Illumination System will change the way you see your patients.

- > **Superior Visualization**—advanced LED light source provides better visualization of the exam area by supplying brighter, whiter light, right where it is needed.
- > **Cordless and Rechargeable**—no cord to get in the way or break during procedures, continuous on time of 80 minutes, recharge time of six hours, lithium ion battery technology. The new cordless illumination system easily fits right into the new KleenSpec® Disposable Vaginal Speculum.
- > **Automatic On/Off**—no more hassle with turning the unit on and off. When the new KleenSpec® Illumination System is inserted into the 590 Vaginal Speculum it is automatically illuminated, and once it is removed the light goes off.
- > **Versatility**—use your new Illumination System as an all-purpose light source as well as with the KleenSpec® Disposable Vaginal Speculum.
- > **Convenient**—no cords to clean, much less heat than traditional lamps.



790 Series Illumination System

Shown with #79920 KleenSpec® Disposable Vaginal Specula Cordless Illuminator with Charging Station and Auxiliary Power Cord



***Welch Allyn is pleased to announce the introduction of our new  
590 Series Vaginal Specula with Cordless Illumination:***



**Small 59000**  
CESS-90634-00



**Medium 59001**  
CESS-90635-00



**Large 59004**  
CESS-90636-00

**590 Series KleenSpec® Vaginal Speculum**

- **Superior Illumination:** newly designed light pipe provides a homogeneous spot and eliminates glare back.
- **Maneuverability:** new wider handle provides improved ergonomics, better balance, and it is much easier to manipulate during exams.
- **Comfortable:** smooth, molded acrylic will not pinch or bind tissue for greater patient comfort during the exam.
- **Compatibility:** the new 590 Series KleenSpec® Vaginal Speculum is compatible with the new cordless illumination system (catalog # 79910) as well as the current illumination system (catalog # 78810).



**790 Series Cordless Illumination System - SKU # CESS-90070-00**

- **Superior Visualization:** advanced LED light source provides improved visualization of the exam area by supplying brighter, whiter light, right where it is needed.
- **Cordless:** no cord to get in the way during procedures.
- **Auto on and off:** the 590 Vaginal Speculum is automatically illuminated on insertion of the illuminator and once it is removed, the light goes off.
- **Versatility:** use your Cordless Illumination System as an all purpose light source as well as with the KleenSpec® Disposable Vaginal Speculum

- **Continuous on time 80 minutes**
- **Fully charge unit in 6 hours**
- **LED life approx. 5500 hours**
- **Battery Lithium Ion**
- **Battery life approx. 6 years**



## New Kleenspec® Vaginal Speculum and Cordless Illumination System Key Features and Benefits

New Kleenspec® Vaginal Speculum	
FEATURES	BENEFITS
590 Vaginal Specula - New profile	Wider handle provides improved ergonomics, better balance and it's easier to manipulate during exams. Scalloped end of light pipe provides a homogeneous spot and eliminates glare back. Will not pinch or bind tissue during exam.
Clear, transparent speculum	Provides excellent visibility of the full exam area.
High quality acrylic	Extremely dependable and reliable. Premium quality. 3X stronger speculum and the yoke is 2X more flexible than low cost substitutes.
Single use disposable	Avoids the DANGER of cross contamination. CONVENIENT – saves time and money associated with cleaning and reprocessing.
Easy one-handed operation	Easy to use operation providing multiple settings for easier exams.
New Kleenspec® Vaginal Speculum Cordless Illumination System	
FEATURES	BENEFITS
LED	Provides brighter, whiter light - 5500 degree Kelvin. Providing true tissue color rendition and enhanced visualization. No heat, no lamp replacement and much more robust.
Superior Illumination	Results in no glare back and a very uniform spot (no shadows or hot spots).
Automatic On/Off	When the new cordless illuminator is inserted into the 590 Vaginal speculum, it automatically illuminates, and when removed, it shuts off.
Cordless & Rechargeable	No cord to get in the way during procedures. Continuous on time of 80 minutes (40 exams) and recharge time 6 hours. Lithium ion battery technology.

# **Frequently Asked Questions**

## **New 790 Series KleenSpec® Vaginal Speculum Illumination System**

What is the continuous on time of a fully charged illumination system?

**Answer:** *80 minutes*

How long does it take to fully charge the illumination system?

**Answer:** *6 hours*

How will I know that the illumination system is fully charged?

**Answer:** *When you put your illuminator in the charger, the indicator on the charger will turn yellow to indicate it is charging. The indicator will return to green to indicate the illuminator is fully charged.*

How many procedures can I perform during the “80 minute on time”?

**Answer:** *Average exam time 2 minutes = 40 exams*

What is the warranty on the battery, charger, and light source?

**Answer:** *1 year*

What is the life of the LED in hours?

**Answer:** *5500 hours*

Am I going to have a warning when the battery charge of the illuminator is low?

**Answer:** *Yes there is a low battery indicator on the illuminator*

While using the aux power cord does the battery charge?

**Answer:** *While the illuminator is lit, the battery does not charge. When you release the switch or remove the illuminator from the speculum, the aux power cord will charge the battery.*

What is the life of the battery?

**Answer:** *3 Years*